

# Cool-TEK Cryo Pneumatic Knee Orthosis w/R.O.M. Hinge

1-877-242-3232 (toll free)  
www.thermotekusa.com

The Cool-TEK Cryo Pneumatic Knee Orthosis with Range of Motion Hinge is designed to provide medial-lateral joint stability, compression and concentrated cold therapy for the knee. The Cool-TEK Knee features dual, adjustable range-of-motion hinges and is medically indicated if the patient has had recent injury to, or a surgical procedure on the knee.

What makes Cool-TEK Knee especially unique is its design. Cool-TEK Knee features plush foam laminate with hook and loop closures, with reusable gel packs encased in a soft fabric lining with nontoxic gel. Cool-TEK Knee combines the benefits of pneumatic compression with cold therapy; featuring a pump bulb that detaches once the desired cold & compression level is achieved. Cool-TEK Knee includes two aluminum ROM hinges with plastic covers and a polycentric flip-top hinge with easy set flexion and extension stops at 0°, 15°, 30°, 60°, and 90°. The Cool-TEK Cryo Pneumatic Knee Orthosis provides excellent medial-lateral control of the knee, and is extremely comfortable for the patient to wear.

ThermoTek's line of Cool-TEK braces are ideal for post-surgical rehab, inflammation, physical therapy, occupational therapy and athletic training.



## SPECIFICATIONS

- Manufactured from plush foam laminate with hook and loop closures
- Fits up to 24" circumference\*
- Billing Code - L1833

\*Measurement taken from leg circumference 4" above mid-patella


## FEATURES & BENEFITS

- Features ROM Hinges with malleable aluminum uprights offering full flexion and extension control with settings at 0°, 15°, 30°, 60° and 90°
- Reusable gel packs feature soft fabric lining and non-toxic gel
- Includes two aluminum ROM hinges with plastic covers and extension straps

## INDICATIONS FOR USE

- ACL repair
- Post total joint surgery
- Meniscal tears
- Patella realignment or during rehabilitation

*The Cool-TEK Cryo Pneumatic Knee Orthosis comes complete, and includes 2 Gel Packs, Valve, Pump & Hose Set*

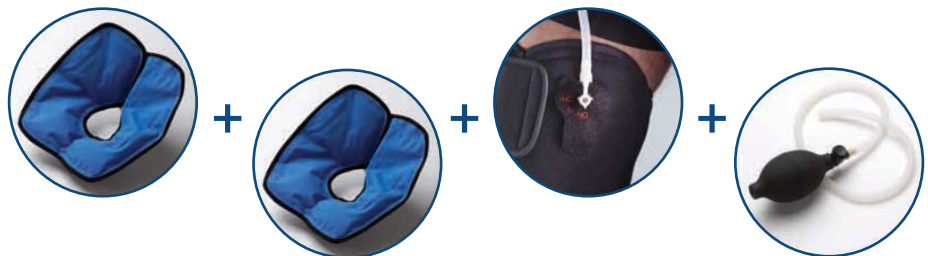

**ThermoTek**<sup>TM</sup>  
1200 Lakeside Parkway, Flower Mound, TX 75028